




## Sentri™ Intersurgical EcoLite™ End Tidal CO<sub>2</sub> Mask

Delivers supplementary oxygen with CO<sub>2</sub> monitoring

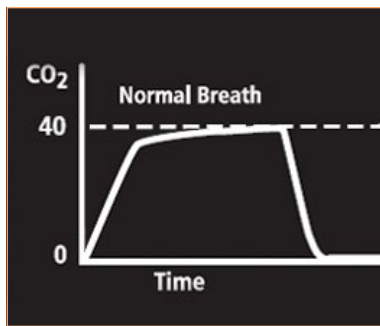

Oxygen and Aerosol Therapy • Variable Oxygen Concentration (Low Flow)



## Intersurgical EcoLite™ Mask

Monitoring ETCO<sub>2</sub> is an invaluable tool for reliable, early warning of respiratory compromise. By using the Senti™ Intersurgical EcoLite™ ETCO<sub>2</sub> mask, no matter where in the hospital setting IV sedation is delivered, both the patient and clinician can take comfort knowing that early detection of inadequate ventilation is made easier with this innovative mask.

The Senti Intersurgical EcoLite Mask is a medium concentration oxygen mask with the addition of a sampling port. This port permits sampling of exhaled carbon dioxide from both the mouth and nose in non-intubated patients during the administration of supplementary oxygen.

### Incurved nose seal

Conforms to different nose shapes, designed to prevent oxygen entering patient's eyes

### No metal nose clip

MRI compatible

### Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

### Soft face seals

Increases patient comfort

### Multi-channel oxygen tube

Oxygen still flows even if tube is kinked

### Luer lock port

Secure ETCO<sub>2</sub> monitoring line connection



### Capnography is vital during sedation



The increased use of conscious sedation has created a need for a device to monitor respiratory depression.

The difference between conscious sedation and general anesthesia is sometimes very small. It is possible during conscious sedation that intravenous sedatives and narcotics administered to alleviate apprehension can result in the loss of consciousness and respiratory obstruction.

The Senti Intersurgical EcoLite Mask is simple to use: attach one end of the ETCO<sub>2</sub> monitoring line to a sidestream capnometer and the other to the luer lock port on the mask.

### A choice of “under chin” or “on chin” positions

Two integral chin seals provide a better fit on a wider range of patient face shapes.

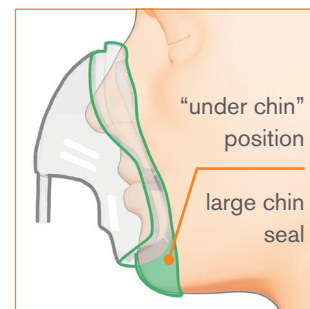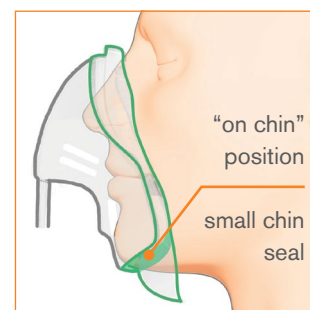

### Comfortable for the environment

Various groups have questioned the use of PVC in medical products and its impact on the environment. Intersurgical's aim is to reduce the environmental impact of its products and processes. The utilization of new materials and the latest manufacturing technology has resulted in the elimination of PVC from the mask shell resulting in a reduced environmental impact.



More about our EcoLite Range



Make an inquiry

[us.intersurgical.com/products/oxygen-and-aerosol-therapy/senti-etco2-medium-concentration-masks](https://us.intersurgical.com/products/oxygen-and-aerosol-therapy/senti-etco2-medium-concentration-masks)

Code	Description	Tube length	Box Qty.
1143015 <i>NEW</i>	Senti Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line, filter and tube	2.1m	30
1142015	Senti Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line and tube	2.1m	30
1141015	Senti Intersurgical EcoLite, adult, mask with oxygen tube	2.1m	30

IS10.24 US • Issue 7 12.18