



ClearLite™ anesthetic face masks

A clear vision for the future 




Anesthesia ▪ Anesthetic Face Masks

ClearLite™ anesthetic face masks

A range of anesthetic face masks which includes seven sizes to meet the clinical needs of all patient types and sizes. The range is available in sizes 0 neonatal to size 6 large adult. Due to the enhanced anatomical design, the ClearLite mask offers superior seal performance with minimum pressure, reducing the risk of unnecessary discomfort or potential trauma.

The superior seal will also assist the clinician in providing optimum manual ventilation of the patient in both the anesthetic and emergency situation.

ClearLite is non PVC and phthalate free which significantly minimizes the impact on the environment. 

Lightweight mask

Provides enhanced grip and comfort in use for the clinician

Anatomical shape

Provides a precise fit for different face types and sizes. Minimum pressure required to achieve an effective seal reducing the risk of inadvertent trauma and nerve damage

Crystal clear body

Provides excellent visibility for confirmation that the patient is well perfused

Reduced mask volume

Reduces potential rebreathing and risks associated with elevated levels of end tidal CO₂

Safe and secure connection

This is achieved between the mask and the breathing system in use. The ISO 22mm tapered connection complies with ISO 5356 Intersurgical airway products

Flexible non inflatable cuff


Provides a comfortable and effective airtight seal









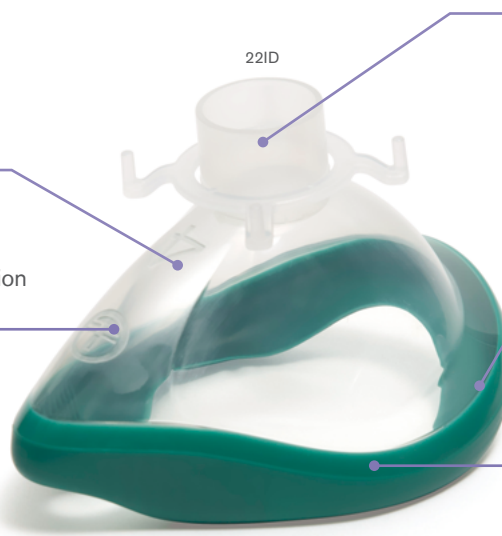





Color-coding

The comfortable, flexible seal is color coded to allow quick and easy identification of the correct sized mask

Low environmental impact

Non PVC and phthalate free



Code	Description	Qty.
 7296000	ClearLite, anesthetic face mask, size 6, extra large adult, red seal, 22ID	20 
 7295000	ClearLite, anesthetic face mask, size 5, large adult, orange seal, 22ID	35 
 7294000	ClearLite, anesthetic face mask, size 4, adult, green seal, 22ID	35 
 7293000	ClearLite, anesthetic face mask, size 3, small adult, yellow seal, 22ID	35 
 7292000	ClearLite, anesthetic face mask, size 2, paediatric, white seal, 22ID	25 
 7291000	ClearLite, anesthetic face mask, size 1, infant, grey seal, no hook ring, 15OD	30 
 7290000	ClearLite, anesthetic face mask, size 0, neonate, blue seal, no hook ring, 15OD	40 

Make an inquiry

Find out more

Watch the video

 Lower environmental impact product

US Lit No. 202.1 • 02.23