




## Air-Guard Clear

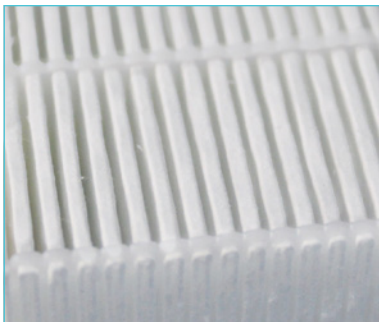

Airway Management ▪ Breathing Filters

## Air-Guard™

### For use in respiratory systems and oxygen concentrators


Intersurgical present the Air-Guard Clear filter, a hydrophobic pleated membrane filter designed for the protection of oxygen concentrator machines and other respiratory equipment.

The Air-Guard Clear filter will function as a microbial filter when used to protect oxygen concentrator machines or when used in breathing systems to protect equipment from contamination.


This product excels as a bacterial and viral filter and can be used as such in a respiratory system providing the maximum period of use is 24 hours, corresponding to the usual correct procedures concerning breathing filters.

**Please note:** this product is contraindicated for use at the patient connection end of the breathing system.

#### Respiratory system protection

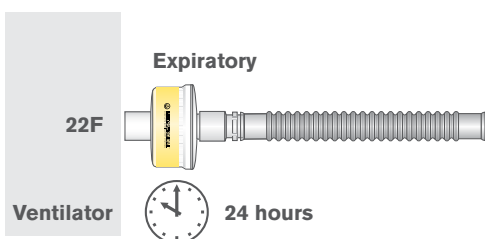

#### Oxygen concentrator protection

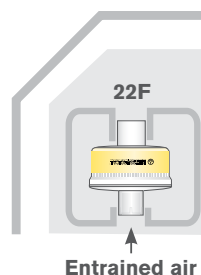

<b>Code</b>	<b>1790000</b>
Box Qty.	50
Filtration efficiency	>99.9999%
Resistance to flow at 30 L/min	0.9 cm H <sub>2</sub> O
Resistance to flow at 60 L/min	1.9 cm H <sub>2</sub> O
Compressible volume	120 ml
Weight	56 g
Connectors	22ID-22OD/15ID
Minimum tidal volume	>360 ml

Make an inquiry

Please think before you print Save energy and paper. If you must print this information sheet please print it double sided.

U.S. Lit No. 82.6 / IS12.14 • Issue 8 07.24



Intersurgical Inc, 2055 Executive Dr. Indianapolis, IN 46241, USA

T: 800-828-9633 F: 315-451-3696 support@intersurgicalinc.com www.intersurgicalinc.com

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • Ireland • Italy • Japan  
Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA