




Cirrus™ 2

A range of single-patient use nebulizers and kits

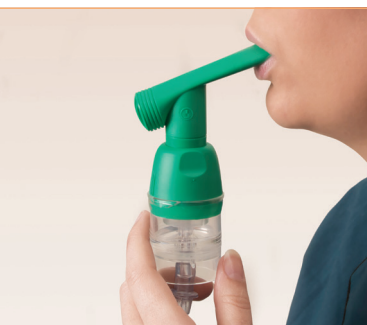

Oxygen and Aerosol Therapy • Nebulizer Therapy




Cirrus™2 nebulizer and nebulizer mask kits

The Cirrus2 nebulizer incorporates a number of features, including a graduated medication cup which provides guidance on how much medication is in the nebulizer. The cup has a quick-release fitting which enables the nebulizer to be attached, as well as detached, with a quarter turn. This speeds up the filling process and provides a more reliable seal between the cup and the nebulizer body.

The Cirrus2 nebulizer is designed to deliver medication for tracheobronchial deposition.

▶ Video available at www.intersurgical.com

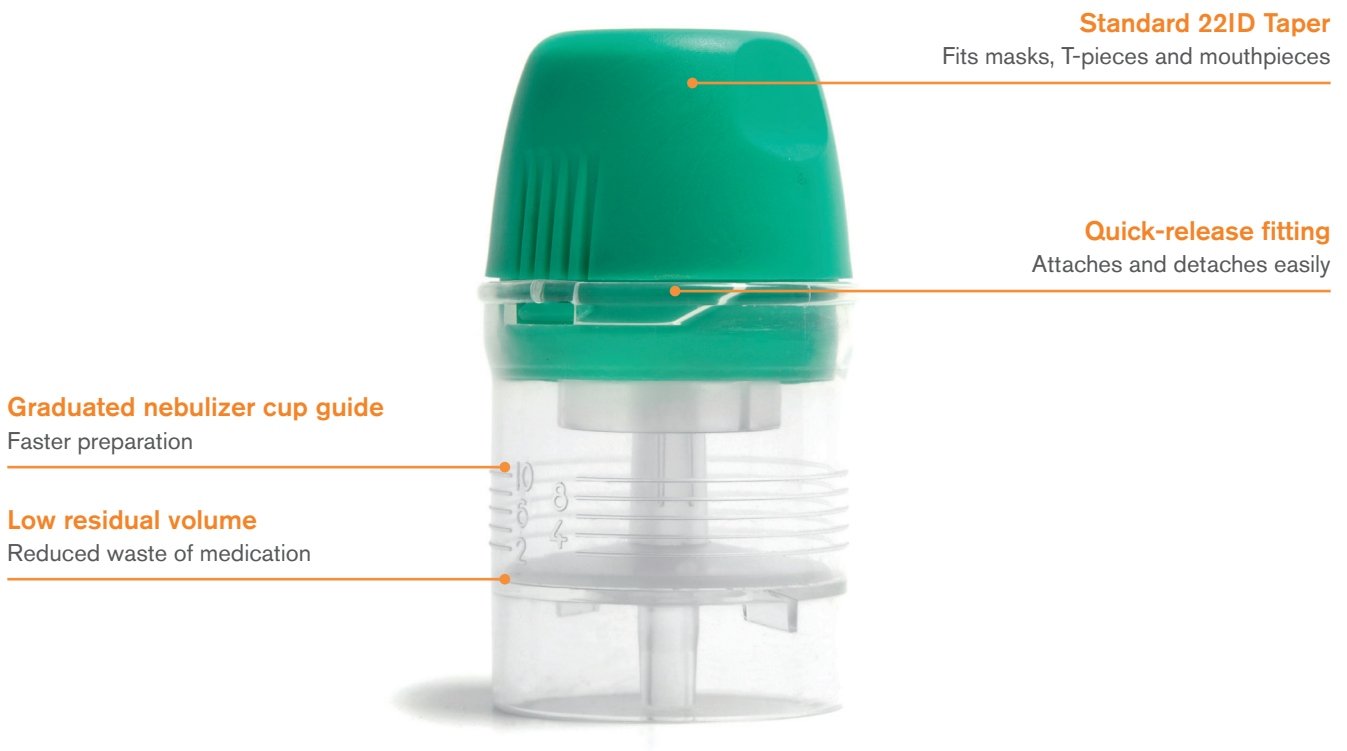



Cirrus™2 nebulizer, Intersurgical EcoLite™ mask kit

Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials into the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, while the second softer thermoplastic elastomer (TPE) material is utilized in the seal, which is in contact with the patients face.

Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilization of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask, reducing our environmental impact.



www.intersurgical.com/info/IntersurgicalEcoLite

Video available at www.intersurgical.com

Incurved nose seal

Conforms to different nose shapes designed to prevent aerosolized medication entering the patient's eyes

No metal nose clip

MRI compatible

Soft face seals

Increased patient comfort


Elastic can be positioned under or over the ears

Below ear position reduces trauma to the top of the ears



A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

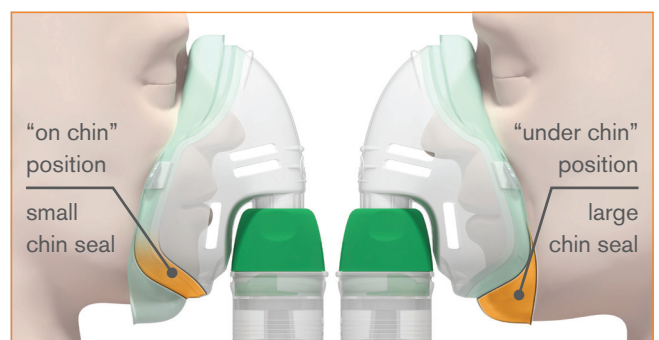



Cirrus™2 nebulizer and nebulizer mask kits

illustrated

Code	Description	Tube length	Qty.
1402004	Cirrus2 nebulizer, single patient use		50
1453013	Cirrus2 nebulizer, adult, Intersurgical EcoLite™ mask kit with tube, single patient use	2.1m	30
1454013	Cirrus2 nebulizer, pediatric, Intersurgical EcoLite™ mask kit with tube, single patient use	2.1m	36
1455000	Cirrus2 nebulizer mouthpiece kit with tube, single patient use	2.1m	35

Nose clip

The nose clip with foam pads are available for use with the Cirrus2 nebulizer.



Code	Description	Qty.
14356000	Nose Clip with foam pads	10

Lower environmental impact product

US Lit No. 266.1 • Issue date 08.22