



DART-Reach™

Laryngo-tracheal mucosal atomization device



Airway Management ▪ Atomization Devices



Quality, innovation and choice

Interact with us




www.intersurgicalinc.com

DART-Reach™ laryngo-tracheal mucosal atomization device

Precise and direct atomization delivery

DART-Reach is indicated for atomizing solutions across the naso, laryngeal and oropharyngeal mucous membranes with an average particle size of approximately 50 µm. Medication can be delivered quickly and effectively for difficult airways with the DART-Reach. The device has a flexible stylet, which allows for a customized approach and can be used in a variety of procedures, including awake intubations, difficult intubations, bronchoscopy and trans-oesophageal echocardiograms.

Our DART-Reach range includes a number of product options, including options with and without a syringe as well as two different lengths of stylet (218 mm and 117 mm). All DART-Reach devices are single patient use and packaged in boxes of 25 units.

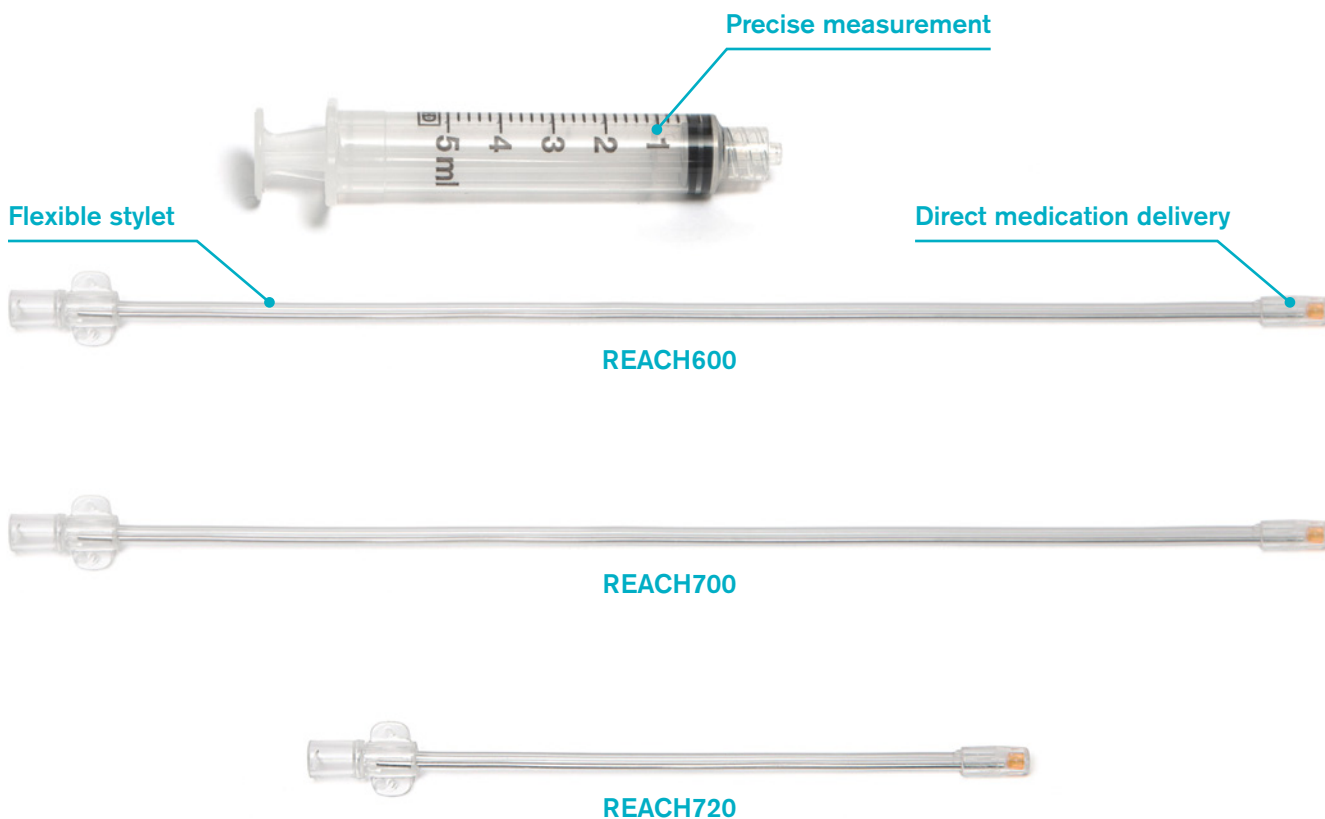

Code	Description	Box Qty.
REACH600	DART-Reach atomization device with 5 ml syringe, 218 mm	25
REACH700	DART-Reach atomization device, 218 mm	25
REACH720	DART-Reach atomization device, 117 mm	25

Make an inquiry

DART-Reach is manufactured by Pulmodyne Inc, and distributed by Intersurgical Ltd

Please think before you print Save energy and paper. If you must print this information sheet please print it double sided.

U.S. Lit No. 309.2 / IS6.23 US • Issue 3 10.24



Intersurgical Inc, 2055 Executive Dr. Indianapolis, IN 46241, USA

T: 800-828-9633 F: 315-451-3696 support@intersurgicalinc.com www.intersurgicalinc.com

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • Ireland • Italy • Japan
Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA